

002-0650-00 Timer Knob Kit for Model M7-001 thru M7-006 Sterilizers

KIT COMPONENTS

Part Number	Description	Qty.
002-0650-00	Timer Knob Kit.....	1
• 053-1008-01	Knob, Timer	1
• 061-0764-00	Label, Timer Knob	1
• 003-1083-00	Installation Instructions.....	1


REQUIRED TOOLS

5/64" Allen Wrench
Alcohol Swab (for cleaning front panel)

INTRODUCTION

The Timer Knob Kit is for replacement of timer knobs on the Model M7-001 thru M7-006 Sterilizers. The kit contains not only the knob but also a self adhesive label that indicates the time span from the OFF position to 30 minutes in 1 minute increments.

INSTALLATION

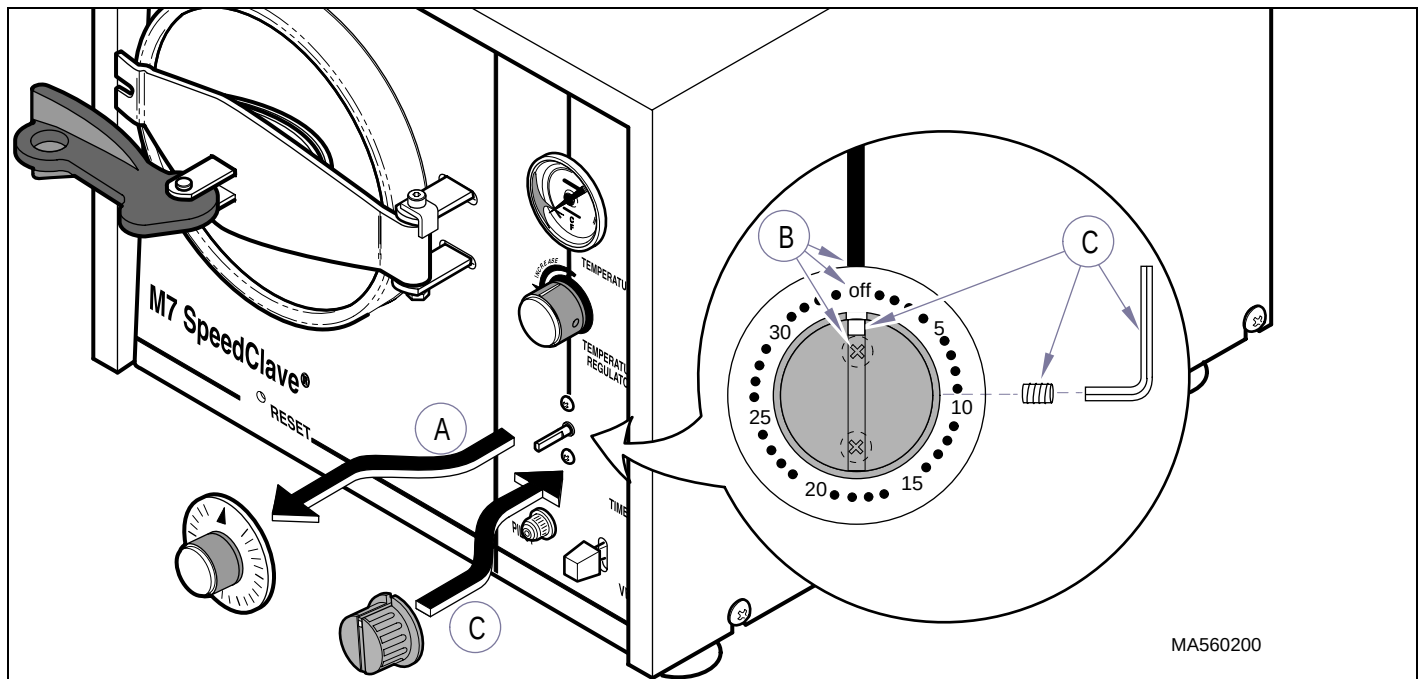

NOTE

Some Units have a pull off timer knob. Other units have an allen screw securing the timer knob to the timer shaft and require a 5/64" allen wrench to loosen the screw.

A. Remove the timer knob.

NOTE

Clean the surface that surrounds the timer shaft before applying the timer label.

B. After removing protective backing from label, position label so the "off" designation is aligned at 12:00, with top mounting screw of the timer or line on front panel then press the label firmly to the surface.

C. With timer shaft in OFF position, align the timer knob's white indicator with the "off" designation on the label. Using a 5/64" allen wrench secure timer knob to shaft and check operation.